

BCF Series Adapters for Leica, Zeiss, Nikon, Olympus Microscopes

BCF series adapters are used to connect C-mount cameras to Leica, Zeiss, Nikon, Olympus Microscopes. The main feature of these adapters is the focus is adjustable, so the images from digital camera and the eyepieces can be synchronous.

BCF series adapters include 0.5× and 0.66× built-in reducing lens respectively. 0.5× adapters are fit for 1/3", 1/2.5", 1/2.3" image sensors. 0.66× adapters are fit for 1/2", 1/1.8" and 2/3" image sensors. Bring you the best field of view for your digital imaging.

Features

1. With the inner focus function, you can synchronize digital images with the eyepiece by just rotating the adapter body slightly.
2. The adapters' head are suit for four major Brands microscope families: Leica, Zeiss, Nikon and Olympus.




BCF-Leica



BCF-Leica0.5x



BCF-Leica0.66x

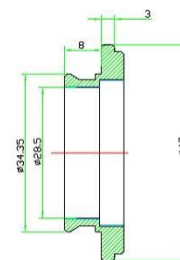

BCF-Leica dimension




BCF-Zeiss




BCF-Zeiss0.5x



BCF-Zeiss0.66x

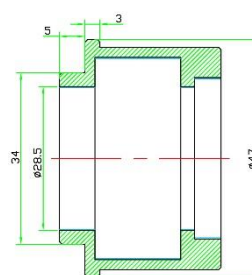

BCF-Zeiss dimension

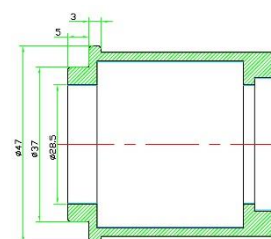


BCF-Nikon

BCF-Nikon0.5x

BCF-Nikon0.66x

BCF-Nikon dimension

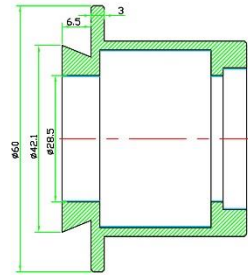

BCF-Olympus

BCF-Olympus0.5x

BCF-Olympus0.66x

BCF-Olympus dimension