


I type Diagnostic Microscope Slides

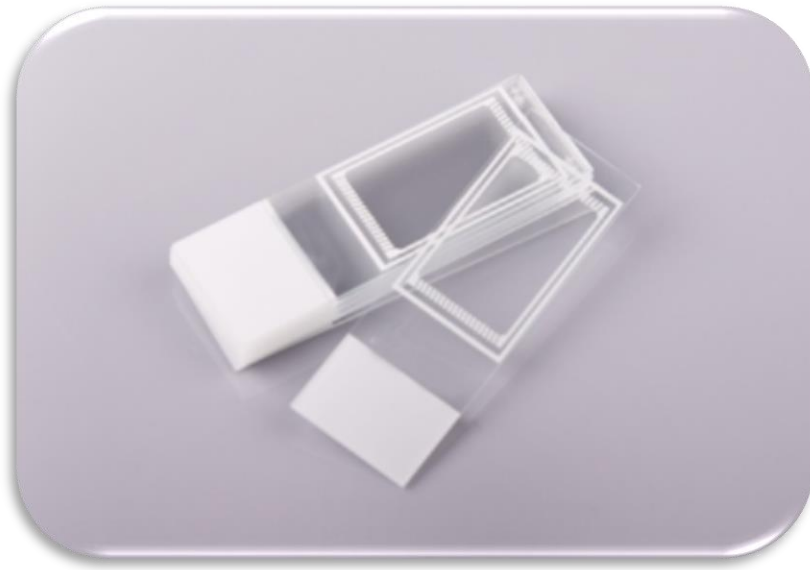

FEATURE

- * The slide surface is coated with PTFE to produce a grid with a thickness of 30-80 μ m. When used, the tissue sections are fixed in the grid, and the antigen repair process of all immune combinations is completed in the grid, which saves a lot of antibodies and reagents.
- * Ideal for manual IHC and automatic IHC with grid slide (such as Biogenex Xmatra Infinity Automatic Staining System).

SPECIFICATION

Item No.	Dimension	Edges	Corner	Packaging	Marking surface	Additional coating
RM7430I	25x75mm 1-1.2mm Thick	Ground Edges	45°	50pcs/box	white	No coating

OPTIONAL

Additional coating: no coating, positive charged coating, hydrophilic adhesion coating.
Customized design and hydrophobic printing color are available.